KOMBILINE

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Content

1. Introduction	5
2. Medical requirements / Applications	6
3. Equipment	6
4. Technical Data	9
4.1. Standard Sizes	9
4.2. Outer measurements of the undercarriage	9
4.3. Air pressure table	9
5. Upholstery / Cover	10
5.1. Attaching and Detaching the back and side upholstery	10
5.2. Removing the cover of the Velcro back and seat upholstery	10
5.3. Detaching the cover of the seat upholstery	11
6. Operation	12
6.1. Adjusting the seat width and depth	12
6.2. Attachment of the seat shell	12
6.3. Removing of the seat shell	13
6.4. Adjustment of the back angle by means of the slide	14
6.5. Adjusting the side supports	14
6.6. Folding the undercarriage	15
6.7. Unfolding the undercarriage	16
6.8. Inclination adjustment of the seat unit	16
6.9. Adjustment of the handle	17
6.9.1. Angle adjustment of the upper part of the handle	17
6.9.2. Height adjustment of the handle unit	18
6.9.3. Inclination adjustment of the handle unit with respect to the vertical position	18
6.10. Adjustments of the middle footplate	19
6.10.1. Vertical inclination adjustment	19
6.10.2. Can be folded backwards	20
6.10.3. Height and depth adjustment	20
6.11. Folding up the anti-tip wheels	20
6.12. Dual rear-wheel drum brakes	21
6.12.1. Applying the parking brake	21
6.12.2. Releasing the parking brake	21
7. Optional Equipment	22
7.1. Safety belt for the pelvis	22
7.2. Belt for the breast	22
7.3. Pressure pad for breast and shoulders	22
7.4. Headrest systems	23
7.5. Arm rests	23
7.6. Abduction pommel (spreading wedge)	24
7.6.1. Abduction pommel integrated	24
7.6.2. Abduction pommel adjustable	24
7.6.3. Abduction pommel swing-type	25
8. Maintenance and care	26
9. Safety instructions	27
10. Declaration of conformity	30
11. Warranty	31
12. Service plan	34

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1. Introduction

Congratulations, you have chosen the combined seating system KOMBILINE by interCo GmbH!

The comfort wheelchair KOMBILINE sets new standards in the care of patients with seating systems. High comfort, simple handling and versatile use are at the forefront of the development. KOMBILINE is used when wheelchair care is medically indicated and low-fatigue sitting over an extended period of time is to be made possible. Patients with paralysis, limb defects or deformities, severe joint contractures or damage should be enabled to participate in social life through own or thirdparty transport. KOMBILINE can be individually adapted to the respective physical limitations and is well suited for reuse thanks to the high degree of standardisation.

KOMBILINE is a two-in-one technical aid. Seat shell and wheelchair are co-ordinated such that a high degree of operating comfort is ensured; simple operating elements contribute towards a clear overview. The aluminium seat shell consists of a seat and back unit joined via physiological pivotal points in the pelvis area. The anatomic basic upholstery ensures an orthopaedically correct sitting position right from the start. KOMBILINE is available in a total of six standardised sizes, which can be individually adapted/additionally equipped to suit the individual requirements of the patient. The choice of technical aid is thus quick, uncomplicated and comfortable.

Should you have additional questions about the seating system, ask us or your specialist dealer. We, the staff of interCo, of course issue the warranty for the procurement of spare parts or new adaptation over years.

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2. Medical requirements /

Applications

Wheelchair care is medically indicated in cases of inability to walk or of severely restricted ability to walk, caused by factors such as:

- paralysis
- loss of limbs
- defects/deformities of the limbs
- severe contractures/damage of the limbs

Differing diagnoses with correspondingly related deficiencies require care with a comfort wheelchair:

- transverse lesion of the cord with paraplegia
- progressive muscular dystrophy
- infantile cerebral palsy
- multiple sclerosis
- hemiplegias
- amputations
- dysmelia/ohokomelia
- neurogenic functional disturbances

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3. Equipment

In its standard version, KOMBILINE consists of a physiologically structured and removable seat and back unit, including seat tilting via gas-pressure springs with large securing levers, parking brakes with locking points, height and angle-adjustable pushing grips, swivelling support rollers, integrated wheelbase extension, tyres 22"/24", drum brake wheels at the rear, tyres 200 x 50 mm at the front, footrest system on the central column with continuous height and angle adjustment via clamp lever, seat and back unit from aluminium with physiological pivotal points in the pelvis area, back-angle adjustment with physiological pivotal points via a gas-pressure spring secured on the back by sliding rail, anatomic seat padding, back and side supports height and width-adjustable, velcrofastening back padding for adjusting the height of the seat, velcro-fastening side padding for regulating the seat width. The system

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is manufactured complete with a cover from highquality material.

Your personal, local specialist advisor will assist you in the planning of an individual solution, taking account of your personal requirements.

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4.2. Outer measurements of the undercarriage

	SIMPLY (22"/24") Seat width 40 cm/45 cm	SIMPLY (12"/16") Seat width 40 cm/45 cm
Length max.	91 cm	86,5 cm
Total width	65 cm/69 cm	57,5 cm/61,5 cm
Height to top of adapter plate	50,5 cm	50,5 cm
Min. heigth push bar (seat angle 90°, push bar 90°)	91 cm	91 cm
Max. heigth push bar (seat angle 90°, push bar 90°)	103 cm	103 cm

Weight (Basic features without seat shell)

Depending on specifications 20 - 22 kg

4.3. Air pressure table

	Max. air pressure (bar)
Front steering wheels	2,5
Back wheels 12"	2,5
Back wheels 16"	3,0
Back wheels 22"	4,0
Back wheels 24"	4,0

Important: The maximum load for SIMPLY is 80 kg!

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5. Upholstery / Cover

The seat shell is equipped with anatomically moulded seat upholstery as well as Velcro back and seat upholstery. The cover is made from high-quality fabric which can be taken off for washing.

5.1. Attaching and Detaching the back and side upholstery

- The back and side upholsteries are attached to the inside of the seat shell with a Velcro fastener. The additional back and side pads are also attached with a Velcro fastener to the inside of the seat upholstery. Bring the pads in the correct position and push to fix them (Pict. 5.1/5.2/5.3).
- To detach the upholstery simply pull them strongly.

5.2. Removing the cover of the Velcro back and seat upholstery

- To remove the cover of the back upholstery, take it out of the seat shell, open the zipper at the lower end and pull the cover off.
- To remove the side upholstery, loosen the Velcro fastener and pull the cover off (Pict. 5.4).

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- To remove the cover of the additional side and back upholsteries, open the zipper and pull the cover off (Pict. 5.5).
- To attach the covers again, proceed in reverse order.

5.3. Detaching the cover of the seat upholstery

In order to detach the cover of the mounted seat upholstery proceed as follows:

- Remove the back and seat unit off the undercarriage.
- Remove the Velcro back and side upholstery as described in Chapter 5.1.
- Loosen the four screws on the bottom side with a hexagon wrench key to disassemble the side upholstery (Pict. 6.1).
- Open the zipper and carefully remove the cover.
- Now you can wash the covers. Please follow the instructions of chapter 9.
- To attach the covers, proceed in reverse order.

Advice: Make sure that all seams are aligned with the edged of the upholstery.

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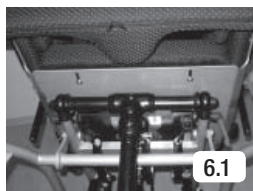
6. Operation

6.1. Adjusting the seat width and depth

The width of the KOMBILINE seat can be enlarged by 5cm by removing the side upholstery (Chapter 5.1.).

The seat depth can also be enlarged by 7 cm.

- Use the four long holes on the bottom of the seat (Pict. 6.1). You will get 3 cm more seat depth. You will get an additional 2 cm of depth by changing the physiological pivotal points and the back mounting suspension to the back drillings (Pict. 6.2).
- By detaching the additional back upholstery (Chapter 5.1) you will get another 2 cm of seat depth.

6.2. Attachment of the seat shell

The seat shell is fixed to the adapter plate of the underframe or wheel chair by means of the trapezoidal adapter (Pict. 6.3). For that purpose perform the subsequent steps and check the firm fit on the adapter plate after the seat shell has been placed.

- First unlock the seat shell fixation underneath the adapter reception by pulling it out (Pict. 6.4)
- Set the rear end of the seat shell adapter to the front of the adapter plate such that the locking is pressed downward (Pict. 6.5).
- Then push the seat shell backwards until it engages audibly (Pict. 6.6).
- Lock the seat shell fixation again.
- Operate the lever of the sliding rail mechanism and lock the back to the push bar until it locks audibly (Pict. 6.7)

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
- Finally, check the seat shell and the seat shell catch for firm fit.


Attention: The placing of the seat shell on the wheel chair may only be performed without patient.

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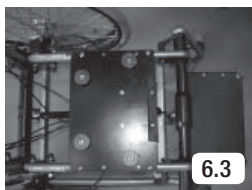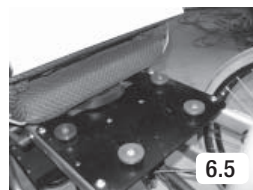6.3. Removing of the seat shell

- First unlock the back rest from the push bar by pulling the latch of the sliding rail mechanism (Pict. 6.8)

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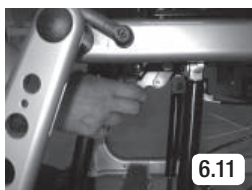
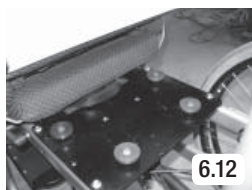
- Pull the lever upwards and the sliding rail will disconnect from the push bar (Pict. 6.9)
- Now loosen the protective device of the seat shell catch (Pict. 6.10).
- Open the seat shell catch by pressing the lever upwards (Pict. 6.11).
- Then you can pull out the seat shell to the front (Pict. 6.12).

Attention: The removal of the seat shell from the wheel chair may only be performed without patient.

6.4. Adjustment of the back angle by means of the slide

Your KOMBILINE is provided with a slide mechanism. Therefore, the adjustment of the back angle conforms to the sliding handle mechanics of the underframe.

- Actuate the hand lever to detach the locking of the backrest mechanics (Pict. 6.12).
- The back angle can be adjusted by pulling and pressing the sliding handle (Pict. 6.13).
- The back angle has been fixed when you release the hand lever.

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6.5. Adjusting the side supports

The side supports are adjustable in height as well as in the lateral dimension to the patient. The adjusting screws are on the backside of the seat shell, two on each support. To adjust please proceed as follows:

- Loosen the adjusting screws with a hexagon wrench until the connection can be moved (Pict. 6.14).
- Now you can adjust the height and the lateral dimensions and fasten the screws again.

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- After adjusting check all connections for firm fit.

6.6. Folding the undercarriage

When folding the undercarriage for transportation, proceed as follows:

- Remove the seat unit (see section 6.3.).
- Loosen the height adjustment of the push bar (see section 6.9.).
- Loosen the screw connections for the handle angle adjustment on the left and right sides of the push bar (see section 6.9.).
- Loosen the catch for the handle folding mechanism by pressing down the central release bar connecting the two locking plates on the lower ends of the two handle upright tubes with your foot while simultaneously swinging the handle forwards.
- The back unit including the handle can now be folded forwards until it rests against the frame.
- Next, remove the rear wheels. In order to do this, lift the frame slightly so that the respective rear wheel is raised off the ground.
- Take hold of the rear wheel spokes and press the button located at the centre of the wheel axle with your thumb.
- Hold the button pressed and remove the wheel from the axle mounting. The wheel brake must be released during this operation.
- Repeat the operation for the second rear wheel.
- If necessary, fold up the anti-tip wheels as described in section 6.11.
- Depending on the version, remove the swivel wheels as described.
- Now, you can load the undercarriage.

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6.7. Unfolding the undercarriage

When unfolding the undercarriage, proceed as follows:

- First of all, check that the parking brakes are released.
- Lift the undercarriage at the rear.
- Take hold of the spokes of one rear wheel and press the button located on the axle with your thumb.
- Hold the button pressed and insert the axle into the axle mounting on the frame.
- Let go of the button and check whether the wheel is fastened securely. It should not be possible to pull it off!

Important: During this check, the button must not be pressed.

- Repeat this operation for the other rear wheel.
- Depending on the version, attach the swivel wheels as described.
- If necessary, turn the anti-tip wheels as described in section 6.11.
- Fold the handle unit back until it engages audibly.
- Next, fasten the loose upper part of the handle in position as described in section 6.8.
- Adjust the desired height of the handle; first lift the handle to the desired height.
- Fasten into position by tightening the quick release connectors as described in section 6.9.
- If clamping is insufficient (the handle can still be moved) proceed as described in section 6.8.
- Attach the seat unit (see section 6.2.).

6.8. Inclination adjustment of the seat unit

The undercarriage (SIMPLY and SIMPLY light) is equipped with a gas spring adjustment for seat unit inclination control. To adjust the seat unit inclination, proceed as follows:

- Apply the parking brakes (see section 6.12.1).
- Pull the gas spring release lever located on the right side of the handle linkage.
- Adjust the inclination of the seat unit by pull-

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ing up or pushing down while keeping the release lever gas pressure spring pulled.

- Let go of the release lever to lock the desired seat inclination position in place.
- If necessary, release the parking brakes (see section 6.12.1)

Important: Never operate the release lever without holding the handle when a patient is seated in the wheelchair.

6.9. Adjustment of the handle

The handle can be adjusted in various ways. This makes an optimal adaptation to the person pushing the wheelchair possible. The following adjustments are possible:

6.9.1 Angle adjustment of the upper part of the handle

To adjust the handle angle, proceed as follows:

- Apply the parking brakes (see section 6.13.)
- Loosen the screw connections of the notched joints on the right and left sides of the handle by turning the clamping lever anticlockwise until the upper part of the handle can be moved.
- Adjust the handle to the desired position.
- Hold the handle in the desired position with one hand and fasten it with the other hand by turning the respective clamping lever until tight.
- Then, fasten the screw connection on the other side by turning the clamping lever until tight.
- When tightening, ensure that the teeth of the notched joints engage correctly.
- If necessary, release the parking brakes (see section 6.12.)

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6.9.2 Height adjustment of the handle unit

A height adjustment of the handle unit is nec-
essary whenever the undercarriage has to be
adapted to the person pushing the wheelchair
or the seat unit inclination has been adjusted.
For height adjustment of the handle, proceed
as follows:

- First of all, loosen the quick release at the
upright tubes by opening the quick release
lever outwards.
- Then, adjust the handle to the desired height
by pulling up or pushing down.
- Fasten in position by tightening the quick re-
lease connectors; fold the lever of the quick
release back into the locked position.
- If the clamping force is insufficient (handle
not fastened securely in place), turn the screw
opposite the quick release lever (with the le-
ver in the open position) clockwise until a
sufficient clamping force is ensured. If neces-
sary, this operation has to be repeated.

6.9.3 Inclination adjustment of the handle unit with respect to the vertical position

In addition to the adjustments described above,
the entire handle can also be tilted backwards by
up to 30°. This may become necessary if, e.g.,
a seat unit with a reclining angle of over 90° is
used. Due to this feature, an extremely small
and even gap between the seat unit back and
the handle is possible. This adjustment is not
aimed at an individual adaptation to the person
pushing the wheelchair. To adjust the inclination
of the handle, please proceed as follows:

- Loosen both clamping levers located in front
of the pivot of the handle folding mechanism
by turning them anti-clockwise.
- Press the coach bolt from the notched guide
by pressing (axially) on the clamping lever
until the inclination adjustment can be freely
adjusted.
- Adjust the handle to the desired position.

Important: For the desired inclination, press
the square section of the bolt into the notch by
pressing on the round head.

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- Check whether the coach bolt is seated tightly in the notch.
- Fasten into position by tightening both clamping levers clockwise.
- Even when folding the handle (e.g. during transportation), the inclination adjustment of the back unit including the handle unit is maintained.

6.10. Adjustments of the middle footplate

When ordered, the middle footplate will already be assembled at the factory. Subsequent individual adjustments are possible. The middle footplate features various adjustment possibilities to guarantee an optimal adaptation to the occupant's needs:

6.10.1. Vertical inclination adjustment

To adjust the middle footplate you have to take the following steps:

- Apply the parking brake (see section 6.12.1.).
- Loosen the middle footplate's screw connections located on the right and left of the upper frame by turning the clamping lever anti-clockwise until the upper part can be moved.
- Bring the middle footplate to the desired position.
- Hold the middle footplate in the desired position with one hand and fasten it with the other hand by engaging the respective clamping lever.
- Then, fasten the screw connection on the other side by turning the clamping lever until tight.
- If necessary, release the parking brakes.

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The footrest is folded up by hand.

6.10.3. Height and depth adjustment

- Loosen the screw at the clamping block of the footrest and adjust the footrest on the middle tube of the middle footplate as required.
- When finished with the adjustments, tighten the screws again.

6.11. Folding up the anti-tip wheels

Two anti-tip wheels are located on the rear of the frame.

Important: In the normal wheeling position, the anti-tip wheels must always point in the direction of travel.

In order to pass over high curbs and during transport (smaller dimensions for transport), it can be necessary to fold the anti-tip wheels upwards. To do so, proceed as follows:

- Apply the parking brakes (see section 6.12.1).
- Stand behind the KOMBILINE.
- First, pull one anti-tip bracket slightly out of the tubular mounting.
- Next, fold the anti-tip bracket in the extended position upwards by 180° so the small wheel points upwards.
- Allow the anti-tip bracket to slide back into the groove in the tubular mounting by means of the spring tension.
- Repeat this operation for the second antitip wheel.
- Release the parking brake as described in section 6.12.1.
- When folding back into original position, proceed in reverse order.

Note: Always make sure that both anti-tip wheels are pointing in the same direction. That is the only way how an effective stabilisation can be ensured! Should you wish to extend the anti-tip wheels, please also make sure that the adjustment is identical on the left and right sides.

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6.12. Dual rear-wheel drum brakes

The undercarriage is equipped with two independently operating drum parking brakes. These act on the two rear wheels. The brake levers are equipped with a three-notch system enabling three different braking levels.

6.12.1 Applying the parking brake

In order to apply the parking brake, proceed as follows:

- Squeeze both brake levers on the handle. Ensure that the brake handle does not return to its initial position (the lock engages audibly); only then, the parking brake is correctly applied.
- Check the braking effect by trying to push the undercarriage.

Important: If the braking action is insufficient, have the brake adjusted by a specialist workshop.

- You can also use the parking brake to make braking easier when going downhill. In this case, pull the brake to the first or second notch position.

6.12.2 Releasing the parking brake

To release the parking brake, proceed as follows:

- Release the lock by operating the notch lever and allow the brake lever to slide into its original position.
- In some cases, tension on the notch lever must be relieved by slightly pulling the brake lever before it is possible to release the brake.

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7. Optional Equipment

7.1. Safety belt for the pelvis

KOMBILINE can be equipped with a belt for the pelvis at a 45°-angle as well as a 70°-angle. A belt is mounted at a 45°-angle after a belt holding device has been mounted to the rear threaded bore and the belt is inserted between backrest and lateral parts. Mounting at the 70°-angle is performed on the lower side of the lateral parts after a belt holding device has been mounted between the two threaded bores and the belt is inserted through the guiding slots provided on the lower side of the lateral parts and moved upwards between lateral parts and seat area (Pict. 7.1).

7.2. Belt for the breast

To mount a belt for the breast KOMBILINE is provided with a mount on the lateral trunk guiding pressure pads (thorax pressure pads). Fasten the belt for the breast to the threaded bores provided.

7.3. Pressure pad for breast and shoulders

KOMBILINE can be equipped with different pressure pads for breast and shoulders. Fastening is identical with the different designs. The lower belts of the pressure pads for breast and shoulders are fastened to the fastening elements of the lateral parts on the rear of the seat shell and moved to the front between backrest and seat area. At the top, the belts are inserted into the tilting buckle on the rear of the seat shell and fixed by the tilting buckle. To fasten the

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patient's belt firmly, the pressure pad for breast and shoulders is loosened at the connecting elements to lay it over the patient. To adapt the pressure pad as to length, use the tilting buckle at the backrest (Pict. 7.2/7.3).

7.4. Headrest systems

KOMBILINE can be equipped with different head rest systems from the series ALULINE (straight or shell shaped). Both are adjustable in height as well as depth. Fix the head rest by loosening the clamping lever of the head rest receptacle on the upper part of the back rest, insert the head rest holder and tighten the lever again. Adjustments are always done by changing the depth of the insertion in the square pipe (Pict. 7.4).

To adjust the depth proceed as follows:

- Turn the fixing lever anticlockwise until the square pipe can be moved (Pict. 7.5).
- Change the square pipe to the desired position.
- Turn the black fixing lever clockwise to fix the square pipe.

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
7.5. Arm rests

KOMBILINE can be equipped with different arm rests from the series ALULINE. Fix the arm rest by loosening the clamping lever of the arm rest receptacle on the right and left side of the seat shell, insert the arm rest holder and tighten the lever again.

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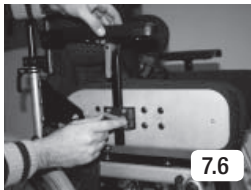Adjustments are always done by changing the depth of the insertion in the square pipe (Pict. 7.6).

To adjust the angle of the arm rest proceed as follows:

- Turn the fixing lever anticlockwise to open the fixture (Pict. 7.7).
- The arm rest angle can be adjusted now.
- Change the square pipe to the desired position and turn the black fixing lever clockwise to fix the square pipe.
- The position of the arm rest will now stay the same.
- Repeat the same procedure on the other side.

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7.6. Abduction pommel (spreading wedge)

Three different abduction pommels are available for KOMBILINE: integrated, adjustable and swing-type.

7.6.1 Abduction pommel integrated

The integrated abduction pommel was affixed during manufacturing of the seat shell. Adjustments can not be made after that. Should adjustments become necessary, please contact your specialist dealer.

7.6.2 Abduction pommel adjustable

The adjustable abduction pommel is affixed at the bottom of the outer side of the hard seat shell. An adjustment is possible. The position of the pommel toward the inner seat can be adjusted. It should be aimed that both knees are centered toward the abduction pommel side.

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To adjust the abduction pommel proceed as follows:

- Loosen the clamp by turning the wing screw anticlockwise. The wing screw is placed under the seat at the mounted holder of the abduction pommel.
- Push the abduction pommel into the correct position.
- Fix the position by turning the wing screw clockwise until the abduction pommel can no longer be moved.

7.6.3. Abduction pommel swing-type

The swing-type pommel is also affixed at the bottom of the outer side of the hard seat shell. For adjustments proceed like with the adjustable abduction pommel (Chapter 7.6.2) To swing the abduction pommel forward proceed as follows:

- Loosen the catch by pressing the red button underneath the abduction pommel (Pict. 7.8).
- Now you can swing back the pommel.
- To bring it back into its original position, move it back until it audibly locks.

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7. Maintenance and care

You can influence the appearance of your valuable home textiles yourself by

- regular cleaning, vacuum-cleaning or careful brushing before dirt gets stuck
- immediately taking up liquids by means of absorbent material – the longer you wait, the more difficult the cleaning becomes! In most cases a mild light-duty detergent is helpful, put it on a cloth (low dose) and use it for dabbing (never rub strongly).

Important: never put this solution directly onto the cover; always use a cloth. Rub off lightly with clean water if necessary.

Please observe: The cloth must be neither chlorinated nor bleached. The tissue should be exposed to direct sunlight as little as possible. All our covering cloths can be washed at 30°. Use a light-duty detergent for washing or gentle chemical purification. Do not dry the cover in a dryer! Should you wish to wash the cover of the seat shell more frequently, a second cover to be used alternately is recommended. Ask your medical supplier for that. Please clean additional parts by means of mild soapsuds and a soft sponge usually used in a household. The basic corpus of the KOMBILINE seat system is provided with a high-quality resisting plastic coating. Clean the latter with soapsuds and a soft sponge. Avoid aggressive or rubbing agents, as they may damage the surface of the plastic coating.

Disinfection

Please use water-based products for the disinfection such as, for example, **Terralin, Quartamon Med** or **Sagrotan original concentrate**. Always observe the manufacturer's instructions for use.

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Service Intervals

SIMPLY LIGHT needs to be checked for functionality and defects by a specialist dealer or us regularly or according to usage, but at least once a year. The necessary checks are listed in the service plan on chapter 12 of this user manual.

9. Safety instructions

- Always comply with the operating and safety instructions of the seating system KOMBILINE.
- A transport of the patient in a car in the seating system KOMBILINE is not permitted at any time! Use special seats meeting applicable requirements for transport in a car. The seating system may only be transported with nobody sitting in it.
- We'd like to point out expressly that persons in seating systems must be secured at any time by means of appropriate holding systems (belts)!
- Always pay attention that the seat shell is firmly connected to the underframe and secured.
- In case of questions related to the KOMBILINE do not hesitate to contact your specialized dealer or InterCo.
- Perform a regular visual inspection of all KOMBILINE connections. If doubts occur, ask your specialized dealer.
- Attachments of accessories or alterations to KOMBILINE performed by others than the manufacturer are not allowed and can result in legal proceedings, if they become known. Further to that, the obligation of warranty by the company of InterCo GmbH expires and the manufacturer can no longer be held liable.
- KOMBILINE is approved for the intended use only.
- Seat units may only be mounted to the underframe in the direction prespecified.

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- Do not use the footboard for entering or leaving the seat system. This should be tilted off or removed.

- To master stairs the patient, the seat system and the underframe must be transported separately. Carrying the underframe including mounted seat system with the patient sitting in it is not allowed! Secure the patient against dropping out before any lifting or lowering operations as well as the adjustment of the sitting inclination.
- To avoid injuries never put your hands between moving parts when performing adjustments.
- Suspending bags is not allowed.
- Never leave the patient without supervision, not even when he/she is secured by belts.
- Should the patient lean out very far or move strongly, the risk of overturning may occur.
- After all adjustment works lock connections which were loosened before again well.

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EU DECLARATION OF CONFORMITY FOR MEDICAL PRODUCTS

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Manufacturer's address: Im Auel 50
53783 Eitorf / Germany

We hereby declare that the product

Name of product: KOMBILINE Wheelchair-Seat-Combination
in 6 standard sizes

Type: Art.-Nr. 37501, 37502, 37503, 37504, 37505, 37506

meets the relevant provisions in the directive specified below:

- Directive 93/42/EWG on medical products, annex I and annex VII

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11. Warranty

In addition to the general terms of business we grant a guarantee of one year for our product KOMBILINE. After returning the guarantee voucher we grant an additional year of manufacturer's guarantee, which means the total guarantee period is two years. Unusable or damaged parts will be repaired within the guarantee period or replaced free of charge. Development-related alterations or wearing parts (e.g. gas pressure damper or covers) are excluded.

We cannot be held liable for damage due to combinations of our product with third parties' products of any kind which may imply considerable risks. We expressly declare that third parties' products approved by us are excluded from that.

Defects resulting from natural wear, excessive strain, damage due to force or with premeditation as well as unintended use are excluded as well.

Subject to technical modifications.

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12. Service manual

Pos.	Area	Check			
		Adjustment / functionality	ok	Damage / deformation	Exchange / replacement
A	Undercarriage				
1	Main frame	Main frame	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Tilting support	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
2	Trapezoidal adaptor		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
3	Steering wheel	Tires	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Tire pressure	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Running performance	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Steering wheel fork	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Suspension on frame	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
4	Driving wheel	Tires	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Tire pressure	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Running performance	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Axle shaft	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Drum brake	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
5	Brakes	Functional test 1	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Functional test 2	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Functional test 3	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Bowden wire	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
6	Adjustment range				
	Gas springs	Damages	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Bolt connection	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
	Seat adjustment	Functional test	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
	Folding mechanism	Folding mechanism	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Arrester	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
	Back angle adjustment	Functional test	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
	Push bar	Height adjustment	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Arrester	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
7	Bolt connections	Size and length	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Completeness	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Damages	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Locktide secure	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Cover flaps in place	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
8	Options	Wheel guard	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Spoke guard	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Arm rests	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		E-Fix	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Middle foot plate	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
9	Optics (Visual inspection)				
	Paint		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
	Cleanliness		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
	Type label		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

Pos.	Area	Check			
		Adjustment / function	OK	Damage / deformations	Exchange / replace
B	Seat shell				
1	Function check	Back adjustment	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Fixing unit	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Arm rest	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Head support	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Abduction wedge	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Therapy table	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
2	Screw connections	Screws tightened	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Rust spots	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
3	Appearance (visual inspection)	Paintwork	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Flock coating	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Cleanliness	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Fabric seams (position / appearance)	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Type plate	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

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Guarantee postcard

If you send us this card, we can assure you an extended guarantee for your KOMBILINE for a total of 2 years. Please fill out this guarantee postcard carefully and completely and return it to us – thank you.

Cost bearer/Medical insurance: _____
Cost bearer's location: _____
Series number: _____
Delivery date: _____
Patient's name: _____
Street: _____
Postcode/Town: _____
Telephone no.*: _____
e-mail*: _____

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